

TRIGLYCERIDES (DST)

(GPO METHOD)

ENZOPAK

Last update 01-2026

Ref.	CC2-TGS.18M, 5x25 ml
	CC2-TGS.18MU, 5x60 ml
	CC2-TGS.18MV, 20x60 ml

INTENDED USE

Reagent kit for quantitative estimation of triglycerides in serum or plasma.

PRODUCT HIGHLIGHTS

- Lipase/GK/GPO-Reagent.
- Very sensitive chromogen.
- Internationally recommended standard giving accuracy of international requirements.

INTRODUCTION

Conventional methods for the estimation of triglycerides have been chemical or enzymatic. In the enzymatic methods, triglycerides are hydrolysed to release glycerol by use of lipase. There are various enzymatic methods to estimate liberated glycerol.

Triglycerides is formulated using Lipo-Protein Lipase (LPL), Glycerokinase (GK), Glycerol-3-Phosphate Oxidase (GPO) and Peroxidase (POD) for quantitative estimation of serum triglycerides. High molar extinction coefficient of the final coloured complex makes the method quite sensitive.


DIAGNOSTIC SIGNIFICANCE

Normally triglycerides, HDL-cholesterol, total cholesterol are estimated, and LDL-cholesterol is calculated. These parameters represent a routine practical aspect of lipid profile which is useful in determination of risk factor or health status of a subject.

Serum triglycerides estimation is an important parameter in the investigation of hyperlipoproteinaemia. Elevated levels may be found in atherosclerosis, diabetes mellitus, glycogen storage diseases like Von Gierke's disease, secondary hyperlipoproteinaemia, alcoholism and nephrotic syndrome.

PRINCIPLE

Lipase hydrolyses triglycerides sequentially to Di & Monoglycerides and finally to glycerol. Glycerol Kinase (GK) using ATP as PO_4 source converts Glycerol liberated to Glycerol-3-Phosphate (G-3-Phosphate). G-3-Phosphate Oxidase (GPO) oxidises, G-3-Phosphate formed to Dihydroxy acetone phosphate and hydrogen peroxide is formed. Peroxidase (POD) uses the hydrogen peroxide formed, to oxidise 4-Aminoantipyrine and DHBS (3,5, dichloro-2-hydroxy benzene sulphate) to a red coloured complex. The absorbance of the coloured complex is measured at 520 nm (500-550 nm) or with GREEN filter which is proportional to Triglyceride concentration.

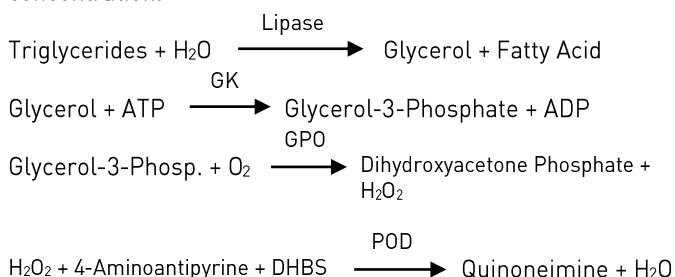

PRESENTATION

All reagents to be stored at 2-8°C

	No. of Bottles		
	5x25ml	5x60ml	20x60ml
• 1 Triglycerides (Enzymes, Chromogen)	5	5	20
• 2 Triglycerides (Buffer)	1	5	20
• Triglycerides Standard (200 mg/dl)	1	1	4

FINAL REAGENT COMPOSITION

ACTIVE INGREDIENTS	Concentration
LPL	≥ 1000 U/L
GK	≥ 800 U/L
GPO	≥ 1500 U/L
POD	≥ 2500 U/L
4 - AAP	0.25 mmol/L
ATP	1 mmol/L
Buffer	100 mmol/L
DHBS	0.5 mmol/L
pH	8.0 ± 0.5 at 25°C
Triglycerides Standard (200 mg/dl)	

Also contains non-reactive fillers and stabilizers.

PRECAUTION

Triglycerides (DST) is for *IN-VITRO* diagnostic use only. Reagent contains Sodium Azide. DO NOT INGEST.

PREPARATION OF WORKING REAGENT

For 5 x 25 ml

Dissolve the contents of one vial of 1 Triglycerides (Powder) with 25 ml of 2 Triglycerides (Buffer). Mix well to dissolve, wait for 5-minutes before use.

For 5 x 60 ml & 20 x 60 ml

Carefully transfer the content of 1 Triglyceride (Powder) into the bottle containing 60 ml of 2 Triglyceride (Buffer). Mix well to dissolve Wait for 5-minutes before use.

REAGENT STORAGE & STABILITY

Triglycerides (DST) reagents are stable at 2-8°C until the expiry date indicated on the label.

Working reagent is stable for 6 months at 2-8°C, when stored in original container protected from light.

SPECIMEN COLLECTION

Fresh, clear fasting serum with no hemolysis should be used. Heparin plasma may be used. No other anticoagulant is suitable. Serum levels are slightly (5mg/dl) higher than plasma levels.

REACTION PARAMETERS

• Type of Reaction	: End Point
• Wavelength	: 520 nm (500-550nm)
• Flowcell Temperature	: 37°C
• Incubation	: 5 min. at 37°C
• Std. Concentration	: 200 mg/dl
• Sample Volume	: 20 μ l (0.020 ml)
• Reagent Volume	: 1.0 ml.
• Zero setting with	: Reagent Blank
• Light Path	: 1.0 cm.

TRIGLYCERIDES (DST)

(GPO METHOD)

ENZOPAK

Last update 01-2026

TEST PROCEDURE

Pipette into Test tubes	Blank	Standard	Test
Working Reagent (ml)	1.0	1.0	1.0
Standard (ml)	-	0.020	-
Sample (ml)	-	-	0.020

Mix and incubate at 37°C for 5 minutes and read absorbance of test and standard against reagent blank at 520 nm (500-550 nm or GREEN filter).

TEST RESULTS

$$\text{Triglycerides (mg/dl)} = \frac{\text{Absorbance of Test}}{\text{Absorbance of Std.}} \times 200$$

To convert (mg/dl) to mmol/lit use the following equation
mmol/lit. = mg/dl x 0.0114

NORMAL VALUES

Male : 65 - 190 mg/dl

Female : 45 - 170 mg/dl

LINEARITY

This method is linear upto 800 mg/dl. For sample values higher than 800 mg/dl, dilute the samples suitably with 0.9% saline and repeat the assay. Apply proper dilution factor to calculate the final results.

REFERENCES

1. FOSSATI P., LORENZO, P., : Serum Triglycerides determined colorimetrically with an enzyme that produces hydrogen peroxide, Clin. Chem 28.2077 - 2080(1982).
2. MCGOWAN, M. W. ARTISS, J. D. STRANBERG, D. R. ZAK, B. A.,: Peroxidase coupled method for the colorimetric determination of serum Triglycerides, Clin. Chem.29, 538-542 (1983)



Regd. Office.3/7, B.I.D.C., Gorwa, Vadodara 390 016 (INDIA)
Web: www.reckondiagnosics.com Ph: +91-265-2281631
Email: mail@reckondiagnosics.com
An ISO 13485:2016 & GMP Certified Company

**RECKON**
DIAGNOSTICS P. LTD.
THE TEST OF CONFIDENCE - SINCE 1989